



Absolute Medical. 1843 W Hubbard St #2A
 Chicago, IL 60622
 Tel: 312 233 2207 Fax: 866 860 9358
 To Order Online: orders@absolutemedical.com
 Our website: absolutemedical.com

Patient Last Name: _____ Patient First Name: _____

Fitter Last Name: _____ Fitter First Name: _____

Date: _____

FarrowWrap® Prescription Form



	Left	Right
Thigh	<input type="checkbox"/>	<input type="checkbox"/>
Leg	<input type="checkbox"/>	<input type="checkbox"/>
Foot	<input type="checkbox"/>	<input type="checkbox"/>
Toes	<input type="checkbox"/>	<input type="checkbox"/>



Mild to Moderate Edema

Moderate to Severe Edema*

*Be cautious when prescribing "Moderate to Severe Edema" for patients with PAD.

If mixed severity, please specify in special instructions below.

	Left	Right
Arm	<input type="checkbox"/>	<input type="checkbox"/>
Hand	<input type="checkbox"/>	<input type="checkbox"/>
Fingers	<input type="checkbox"/>	<input type="checkbox"/>



Special Instructions: _____

Related Diagnoses: _____


Patient Last Name: _____ Patient First Name: _____

Fitter Last Name: _____ Fitter First Name: _____

Date: _____

FarrowWrap® Thigh High OTS and Custom Order Form

Note: All measurements in cm

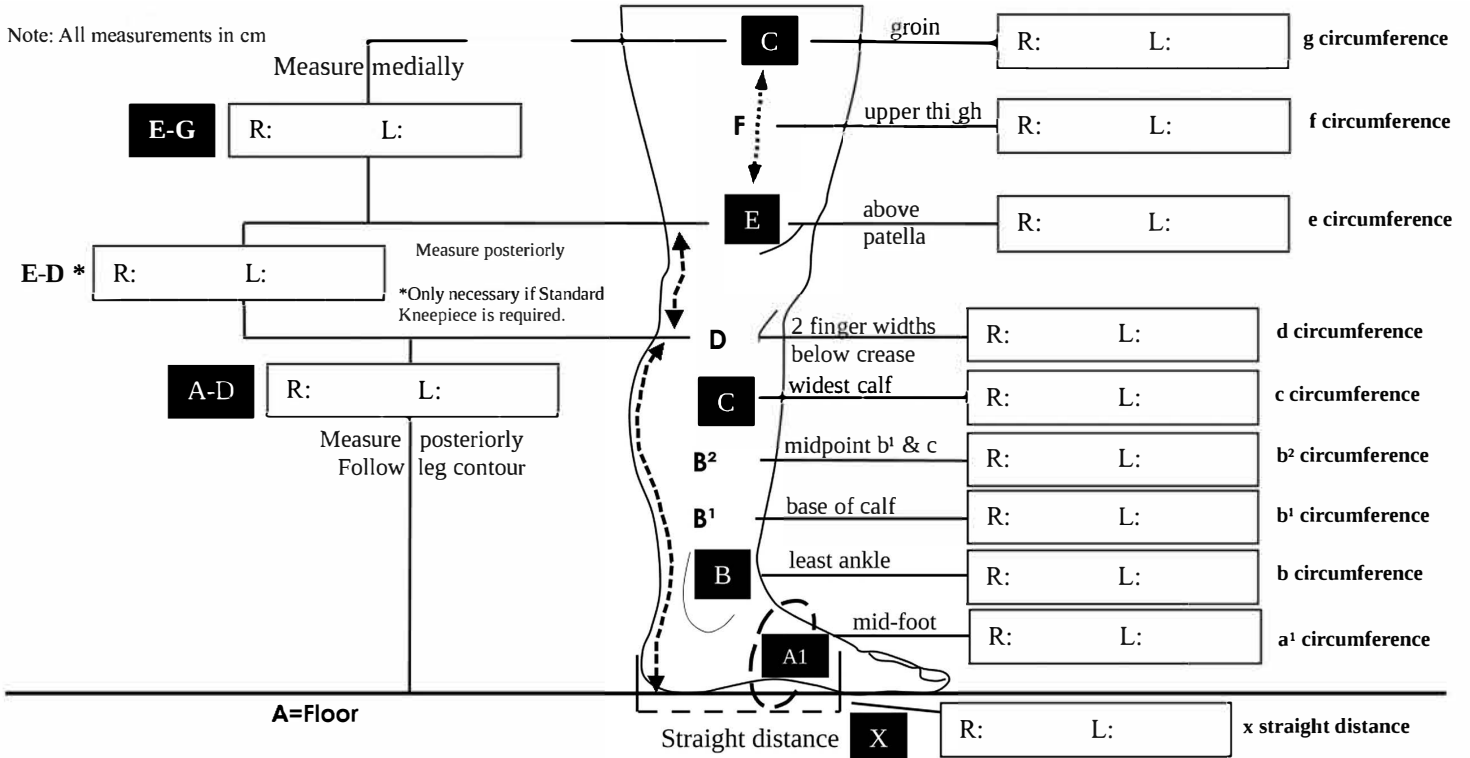

Please place a numeric value in the white box to indicate items needed.	Off-The-Shelf Thighpiece Size Chart (Only comes with wrap-around kneepiece)									
	XSmall	Qty	Small	Qty	Medium	Qty	Large	Qty	XLarge	Qty
(G) Groin	60–70 cm		65–75 cm		70–80 cm		75–90 cm		85–100 cm	
(E) Above Knee	47–52 cm		53–59 cm		60–66 cm		67–73 cm		74–80 cm	
(E-G) Short FW_-O-TS	16 – 20 cm		16 – 20 cm		16 – 20 cm		16 – 20 cm		16 – 20 cm	
(E-G) Reg FW_-O-TR	21 – 26 cm		21 – 26 cm		21 – 26 cm		21 – 26 cm		21 – 26 cm	
(E-G) Tall FW_-O-TT	27 – 31 cm		27 – 31 cm		27 – 31 cm		27 – 31 cm		27 – 31 cm	
Please place a numeric value in the white box to indicate items needed.	Off-The-Shelf Legpiece Size Chart									
	XSmall	Qty	Small	Qty	Medium	Qty	Large	Qty	XLarge	Qty
(C) Widest Calf	36 – 43 cm		42 – 50 cm		48 – 58 cm		53 – 63 cm		58 – 68 cm	
(B) Ankle	21 – 25 cm		25 – 30 cm		30 – 36 cm		36 – 42 cm		42 – 50 cm	
(A-D) Reg FW_-O-LR	33 – 37 cm		35 – 39 cm		37 – 41 cm		39 – 43 cm		39 – 43 cm	
(A-D) Tall FW_-O-LT	38 – 41 cm		40 – 43 cm		42 – 45 cm		44 – 47 cm		44 – 47 cm	
Please place a numeric value in the white box to indicate items needed.	Off-The-Shelf Footpiece Size Chart									
	XSmall	Qty	Small	Qty	Medium	Qty	Large	Qty	XLarge	Qty
(A¹) Mid-Foot	22 – 24 cm		25 – 27 cm		28 – 30 cm		31 – 34 cm		35 – 40 cm	
(X) Reg FW_-O-FR	16 – 17 cm		17 – 18 cm		19 – 20 cm		20 – 21 cm		22 – 23 cm	
(X) Long FW_-O-FL	18 – 19 cm		19 – 20 cm		21 – 22 cm		22 – 23 cm		24 – 25 cm	

Options	
FABRIC	✓
Classic	
LITE	
STRONG	
When a footpiece, legpiece and thighpiece are purchased together, 1 pair of Terry Cloth liners are included automatically free. Other liners may cost extra. See liner form & pricing for details.	

Custom Component	Thighpiece	<input type="checkbox"/> Wrap-around <input type="checkbox"/> Standard	Legpiece	FW_-C-L	Footpiece	FW_-C-F	Extend footpiece length to base of toes
Quantity	Right: <input type="text"/>	Left: <input type="text"/>	Right: <input type="text"/>	Left: <input type="text"/>	Right: <input type="text"/>	Left: <input type="text"/>	L R BL: <input type="text"/>

Other: _____

Directions for Thigh High Garments

Please note that OTS sizes can accommodate about a 20% reduction in edema, but only a 5% increase.

Fill in the Patient / Billing & Shipping Information areas **completely**. Next, choose the FarrowWrap® version; Classic, LITE or STRONG. Classic and STRONG are generally reserved for more severe edema cases. LITE is more appropriate for mild to moderate edema. On the measurement page, part of the sku number is represented by a “_” which represents either a “CL” for Classic, “LT” for LITE or “ST” for STRONG. The appropriate SKU number can be crossreferenced with the Pricing Schedule to calculate the price.

Only circumferences at points **A**, **B**, **C**, **E**, & **G** and length measurements **X**, **A-D**, & **E-G** are necessary to determine if the patient can fit into an **OTS** (Off-The-Shelf or ready-made) garment. Length measurements for the legpiece are measured along the posterior contour of the leg and length measurements for the thighpiece are measured medially and following contours around any lobes. If any lobes are present, drawings and/or pictures (with the patient’s permission) will help. Follow these instructions on how to properly measure.

FarrowWrap® Thighpiece

1. Obtain the following 3 measurements:
 - a. Upper thigh circumference (at **G**) just below the groin.
 - b. Distal thigh circumference (at **E**) right above the knee.
 - c. Length of thigh from just above the knee joint to the groin. Measure medially. **E-G**.
2. Now, match the upper and distal thigh circumferences to a respective size. If the upper thigh AND distal thigh circumferences do not both match a single size, or if the patient has a large lobule, then a custom thighpiece is necessary along with all measurements.
3. If the circumferences do match an OTS size, then compare **E-G** to the OTS Short, Regular and Tall length measurements. If the **E-G** length falls within the specified ranges, then in the space next to the appropriate OTS length enter the number of thighpieces you would like in that given size. If the measured length is not within the ranges, a custom thighpiece will be necessary. Enter the number of custom thighpieces required at the bottom of the page in the **Custom Component** section.
4. All OTS thighpieces come with a wrap-around kneepiece. Custom thighpieces come with a standard kneepiece or a wrap-around kneepiece. See our catalog for descriptions. If choosing a custom thighpiece, enter the number of pieces needed for the RIGHT and/or LEFT along with the kneepiece desired, Wrap-Around or Standard.
5. We recommend biker’s shorts or GarmentGrip™ to help hold thigh pieces in place on problematic patients.

FarrowWrap® Legpiece and Footpiece

1. See FarrowWrap® A-D (Knee High) Order Form Instructions

Note

The measurements used by this form allow for at least a 20% additional reduction in leg and thigh circumferences (edema reduction), but can only accommodate about a 5% increase. There is room for greater variance in calf circumferences than for ankle circumferences. Patients may fit into an OTS garment even though this form may indicate otherwise, but we cannot guarantee it.

Liner Options

All footpiece/legpiece/thighpiece combinations will be issued a free pair of Terry Cloth AG Liners automatically. To order additional liners, please use the “Liners Order From” section. We recommend at least 2 pairs of liners per extremity. Terry Cloth Liners come in 3 circumferences, S <40cm, M 40-70 cm, & L 70-120 cm.

Other Options

Hook Stays can be used to help prevent FarrowWrap® legpieces from sliding on highly triangular legs and can be easily attached to the garment. Stays are available in small (16 cm length), medium (21 cm length), and large (27 cm length).