



lymphology measuring for custom-made flat-knit lower extremity

effective October 1, 2017



# Notes on taking measurements for custom-made flat-knit stockings


- Flat knit compression garments should not be measured until the best possible decongestion effort has been achieved. Circumference (¿) and length (¿) measurements are taken in a distal to proximal sequence (starting at the toes and working up the leg towards the waist).
- It is essential to mark measuring points on the leg so that the circumference and length measurements are taken at the same point.
- The amount of tension applied must be determined by a person experienced in bandaging techniques and depends on edema and tissue characteristics, sensitivity to pain, fibrosis, location of the edema, mobility, and strength of the patient.
- Generally flat-knitted garments are measured under tension with the exception of cY. According to the measuring point we distinguish:

Skin measurements = loose; Measuring to the first resistance = measuring under tension without force; Full Tension = absolutely compression of the tissue.

# **Knee high**

All length measurements are contoured along the medial (inner) side of the leg

# Foot measurement points (Y, A, A1, Z)

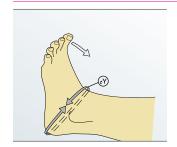

# Circumference (c) Y

Measurements taken over the heel and arch in dorsal flexion with tendons taut. Take physical dimensions without tension. Take soft tissue projections into account.

# Length (ℓ) A


Distance from heel to base of little toe (length of outer side of foot), including the joint.

# Length (ℓ) A1

Distance from heel to base of big toe (medial length), including the joint.

# Length (ℓ) Z

Closed toe: distance from heel to end of longest toe. (Not needed for open toe garments, except open toe seamed toe caps.  $\ell$ Z needed for ALL seamed toe cap garments.)

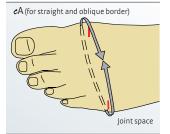

## Circumference (c) A

Apply tape measure around the metacarpophalangeal joints (first joints of toe).

# Measuring point B: just above the malleolus (ankle bone)




# Length (ℓ) B

From sole of foot to above (proximal to) the ankle .

#### Circumference (c) B

Measure circumference, depending on indication, tightly if necessary.

## Measuring point B1: at the base of the calf

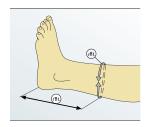


#### Length (ℓ) B1

From sole of foot to base of the calf.

# Circumference (c) B1

Measure circumference, depending on indication, tightly if necessary.

# Measuring point C: at the largest calf circumference

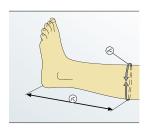


# Length (ℓ) C

From sole of foot to largest calf circumference.

#### Circumference (c) C

Measure circumference, depending on indication, tightly if necessary.

# Measuring point D: at the head of the fibula

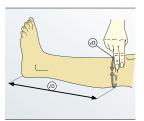

# Length (ℓ) D

From the sole of the foot to the head of the fibula (approx. two finger widths below kneecap).


# Circumference (c) D

Measure circumference, depending on indication, tightly if necessary. Exception: AD stocking = Body dimension.


# Thigh high

All length measurements are contoured along the medial (inner) side of the leg except (I)K2 which is contoured on the posterior (back) side of the leg.

## Measuring point E: at the middle of the kneecap

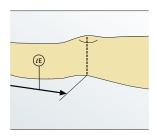


Length ( $\ell$ ) E From sole of foot to middle of kneecap.

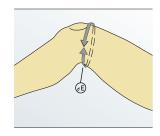

# Circumference (c) E

In general, measurements are not made tightly when the leg is bent (angle 30° - 45°), take circumference measurement one finger width proximal (above) to the back of the knee, and two finger widths proximal (above) to the knee cap. Do not place tape in the back of the knee or folds of skin.

#### Measuring point F: at the middle of the thigh

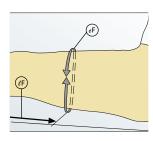

Length (*l*) F
From sole of foot to middle of thigh.


Circumference (e) F Depending on the indication, the circumference is generally measured fairly tight.





Length (ℓ)G

AG stocking dimension: sole of the foot to largest point (inner aspect) below the transverse buttock fold (including top band).

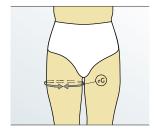

Circumference (c) G

Depending on the indication, the circumference is generally measured fairly tight at the widest part of upper thigh (not always the top of the thigh.)

### Measuring point K1: at the pubic bone



Length ( $\ell$ ) K1 From sole of foot to the pubic bone.

# Measuring point K2: at the gluteal fold



Length ( $\ell$ ) K2 From sole of foot to gluteal fold.


Circumference (e) K Measure around both legs at the gluteal fold.

Obligatory measurement for medi!

# Waist high

All length measurements are contoured along the body except "H" and "T" which are straight lengths.

# Measuring point H: widest point of the pelvis

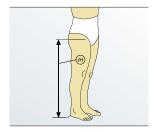


**Length (***t***) H**From sole of foot to widest point of pelvis.



**Circumference (***e***) H**Circumference of pelvis at the widest point, or buttocks, measure tightly depending on indication.

# Measuring point T: at the waist

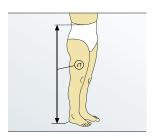


Length (/) T From the sole of the foot to the top border of the garment (patient-specific).




Circumference (c) T Circumference of torso, measured under moderate tension.

### Measuring points for front (¿) TK1 and back (¿) TK2 top border of the garment

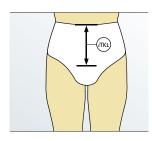

Length (?) TK1
Front length TK1: measure vertically next to body, at the level of the pubic bone to the top border of the garment.

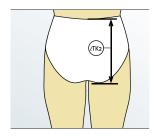

Length (4) TK2

Back length TK2 measure

vertically near the body, at the level
of the transverse gluteal fold to the
border of the garment.

# Measuring a suspensory

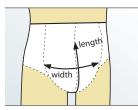


Fig. A

# Length & width

Measure the width and length of the suspensory, as marked in Figures A and B - this measurement can be taken over the underwear.

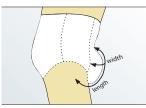



Fig. B




Absolute Medical. 1843 W Hubbard St #2A Chicago, IL 60622 Tel: 312 233 2207 Fax: 866 860 9358 To Order Online: orders@absolutemedical.com Our website: absolutemedical.com

Patient Last Name:	Patient First Name:
Fitter Last Name:	Fitter First Name:
Date:	

# LOWER EXTREMITY MEASURING & ORDER FORM: CUSTOM-MADE FLAT-KNIT

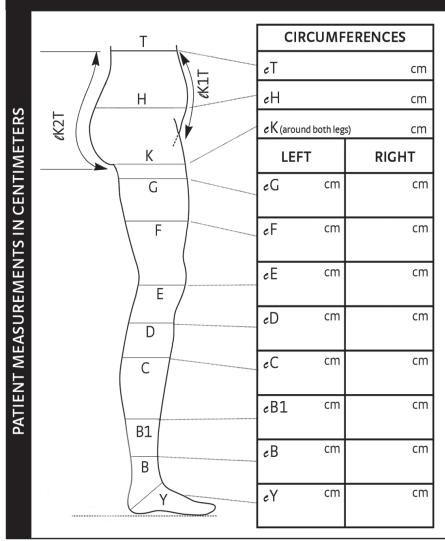

LENGTHS (landmarks to floor)					
ℓT (straight	line)	ℓK1T			
ℓH (straigh	nt line)	ℓK2T			
LEFT		RIGHT			
€G	cm	cm			
<i>ℓ</i> F	cm	cm			
€E	cm	cm			
ℓD	cm	cm			
eС	cm	cm			
<i>€</i> B1	cm	cm			
ℓB	cm	cm			

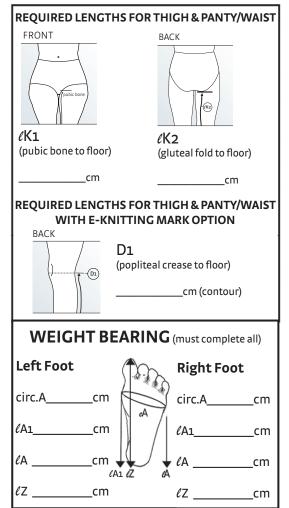



Absolute Medical. 1843 W Hubbard St #2A Chicago, IL 60622
Tel: 312 233 2207 Fax: 866 860 9358
To Order Online: orders@absolutemedical.com
Our website: absolutemedical.com

atient Name:	 Date:	
itter Name:		

# LOWER EXTREMITY MEASURING & ORDER FORM: CUSTOM-MADE FLAT-KNIT

		emaii: orders@absoluter						
ATERIAL	☐ mediven® 550 ☐ mediven® mondi ☐ mediven® sensoo (cclii)		OPTIONS	FOOT  Open toe Closed toe:  Tricot (half compression) *default		OBLIQUE ENDING  ☐ Standard ☐ Straight	☐ Steep	
STYLE/MATERIAL	NOTE: ONLY ONE STYL  ☐ Knee high ☐ Waist high ☐ One leg Waist high* ☐ Capri shorts (T to C)* ☐ Men's leotard	☐ Thigh high☐ Thigh w/waist attachmer		MONDI SILVER  Y-C  Y-D  Y-G	Gluteal shaper (tricot gusset only) G-D G-G Hallux ease (closed toe or seamless toe cap only)  NITTING OPTIONS  Y knitting mark at ankle E knitting mark at knee  LICONE TOPBAND None		□ ADDED LENGTH IN FRONT (include \( \epsilon \text{1} \) and \( \epsilon \text{2} \) measurements below) Needed for all thigh high stockings □ GUSSET: □ SUSPENSORY:	
COMPRESSION	□ with fly  PANTY SECTION  CCLI: □ 18-21mmHg  CCLIII: □ 34-46mmHg	(slight compression panel)  CCLII: □ 23-32mmHg  CCLIV: □ 49-60mmHg **		SILICONE TOPB  None			Tricot (standard Netting Compressive	) Tricot Netting Compressive Width cm
COMPF	LEFT LEG CCLI: ☐ 18-21mmHg CCLIII: ☐ 34-46mmHg	CCLII: ☐ 23-32mmHg CCLIV: ☐ 49-60mmHg**		☐ Profile (honeycomb pattern, 5cm wide only) ☐ Beaded: ☐ Small (2.5cm wide) ☐ Wide (5cm wide) ☐ Peony (beaded, 5cm wide only)		Length cm  LYMPHPADS*  Removeable (standard)  Width (cm) Length (cm)	rd)	
_	RIGHT LEG  CCLI: 18-21mmHg  CCLIII: 34-46mmHg	CCLIV: ☐ 23-32mmHg CCLIV: ☐ 49-60mmHg**		POSITION Lateral/along oblique border		ANTI-SLIP SILICONE DOTS  □ 10 X 4.5 cm	OR ☐ Silk Lining Width(cm) Length(cn	Length(cm)
OR QTY	Left  Standard Colors  Caramel Black	Right Pair  Trend Colors*  Anthracite Magenta		Anterior/vertically over "E"  Posterior/rear over seam	□ 8 X 5 cm □ 8 X 5 cm	□ 10 X 4.5 cm □ 10 X 4.5 cm		or   Medial   Lateral
COLOR	☐ Cashmere* ☐ Sand	☐ Aqua ☐ Mint-Gree ☐ Blue-Jeans ☐ Moss-Gree ☐ Cherry-Red ☐ Navy ☐ Coral		Sole of the foot  WAISTBAND		□ 6 X 4.5 cm	SWAROVSKI CRY	B to D B to G
DESIGN	mediven 550 Fashion Elements* Must pick1color & 1 pattern Colors Patterns  Berry Animal Brown Crosses Gray Ornaments  mediven 550 Design Elements  Stars Ribs Pyramids			□ Adjustable (Standard) □ Velcro □ Silicone beaded (5cm wide) □ Knitted ending  LEVAMED SILICONE ANKLE PADS*  Add 1 cm to cY measurement per pad.  Left: □ Inner □ Outer □ Permanent □ Removeable Right: □ Inner □ Outer □ Permanent □ Removeable		Pattern:  Beauty Brilliance Romance		
	*Not available in mediven s **mediven 550 only NOTES:	ensoo. Trend colors and fashion	elements	require an extra 5 d	ays for produc	tion.		