




lymphology
measuring for custom-made
flat-knit toe caps

effective September 1, 2017

Take the following measurements while patient is weight bearing.

Measuring point A-A1

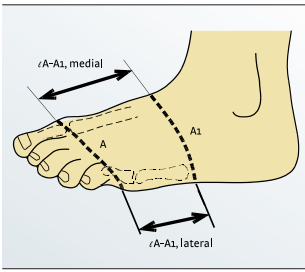


Length (ℓ) A

Medial length: from base joint of big toe to start of arch.

Length (ℓ) A1

Lateral length: From base joint of little toe to start of arch.

Measuring point A & A1

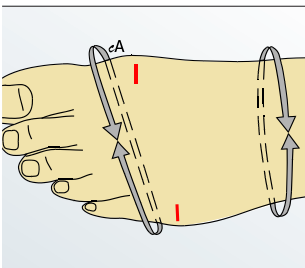


Circumference (c) A

Measure straight across base joint of little toe and base joint of big toe .

Circumference (c) A1

Measure straight across middle of foot (final length).


Measuring point toes 1-5

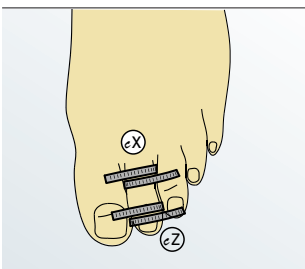

Circumference (c) X and Z

Generally measure dimensions X and Z to the millimeter and without tension. The X measurement should not be smaller than the Z measurement, and, where necessary, is adjusted to the Z measurement. Ok to add up to 5mm to cX measurement.



Measuring point Toes 1-5

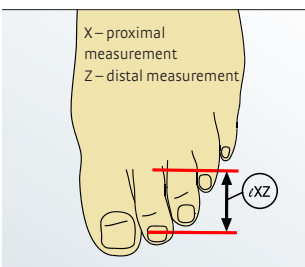

Length (ℓ) XZ

Shortest length from web X to desired length Z.

mediven® Flat-knit Toecaps

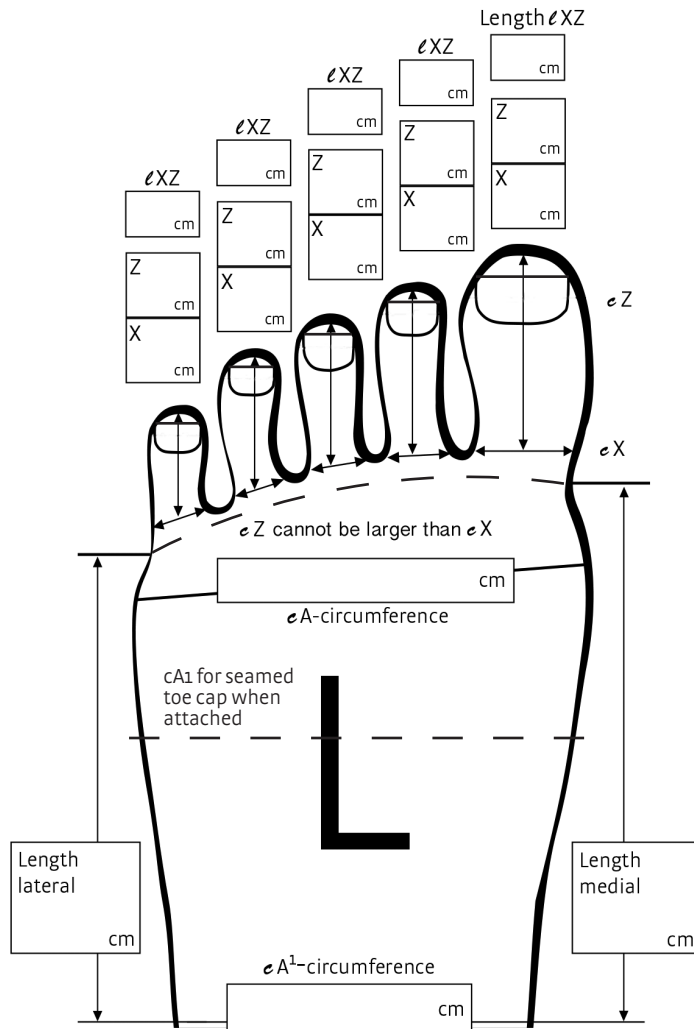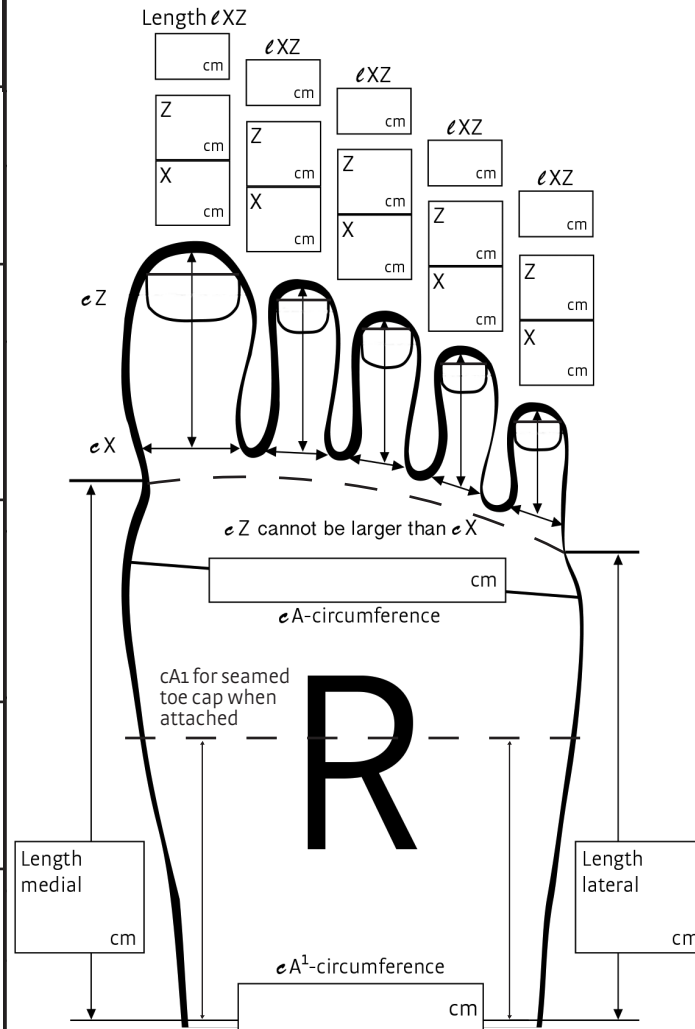	<table border="1" style="width: 100%; border-collapse: collapse;"> <tr> <td style="width: 50%;">Foot</td> <td style="width: 50%;">Material</td> </tr> <tr> <td>_____ Left quantity</td> <td><input type="checkbox"/> 550</td> </tr> <tr> <td>_____ Right quantity</td> <td><input type="checkbox"/> mondi</td> </tr> </table>	Foot	Material	_____ Left quantity	<input type="checkbox"/> 550	_____ Right quantity	<input type="checkbox"/> mondi	
Foot	Material							
_____ Left quantity	<input type="checkbox"/> 550							
_____ Right quantity	<input type="checkbox"/> mondi							
<p>Style (Choose One)</p> <p><input type="checkbox"/> Seamless toe cap (as attachment) must be same compression as stocking</p> <p><input type="checkbox"/> Seamed toe cap (as attachment) can be a different compression than stocking.</p> <p><input type="checkbox"/> Individual toe cap</p>								
<table border="1" style="width: 100%; border-collapse: collapse;"> <tr> <td style="width: 50%;">Standard Colors</td> <td style="width: 50%;">Trend Colors*</td> </tr> <tr> <td> <input type="checkbox"/> Caramel <input type="checkbox"/> Cashmere <input type="checkbox"/> Black <input type="checkbox"/> Sand </td> <td> <input type="checkbox"/> Anthracite <input type="checkbox"/> Aqua <input type="checkbox"/> Blue-Jeans <input type="checkbox"/> Cherry-Red <input type="checkbox"/> Coral <input type="checkbox"/> Magenta <input type="checkbox"/> Mint-Green <input type="checkbox"/> Moss-Green <input type="checkbox"/> Navy </td> </tr> </table>			Standard Colors	Trend Colors*	<input type="checkbox"/> Caramel <input type="checkbox"/> Cashmere <input type="checkbox"/> Black <input type="checkbox"/> Sand	<input type="checkbox"/> Anthracite <input type="checkbox"/> Aqua <input type="checkbox"/> Blue-Jeans <input type="checkbox"/> Cherry-Red <input type="checkbox"/> Coral <input type="checkbox"/> Magenta <input type="checkbox"/> Mint-Green <input type="checkbox"/> Moss-Green <input type="checkbox"/> Navy		
Standard Colors	Trend Colors*							
<input type="checkbox"/> Caramel <input type="checkbox"/> Cashmere <input type="checkbox"/> Black <input type="checkbox"/> Sand	<input type="checkbox"/> Anthracite <input type="checkbox"/> Aqua <input type="checkbox"/> Blue-Jeans <input type="checkbox"/> Cherry-Red <input type="checkbox"/> Coral <input type="checkbox"/> Magenta <input type="checkbox"/> Mint-Green <input type="checkbox"/> Moss-Green <input type="checkbox"/> Navy							
<table border="1" style="width: 100%; border-collapse: collapse;"> <tr> <td style="width: 50%;"> mediven 550 Fashion Elements* <i>Must pick 1 color & 1 pattern</i> Colors <input type="checkbox"/> Berry <input type="checkbox"/> Brown <input type="checkbox"/> Gray </td> <td style="width: 50%;"> mediven 550 Design Elements <input type="checkbox"/> Animal <input type="checkbox"/> Crosses <input type="checkbox"/> Ornaments <input type="checkbox"/> Stars <input type="checkbox"/> Ribs <input type="checkbox"/> Pyramids </td> </tr> </table>			mediven 550 Fashion Elements* <i>Must pick 1 color & 1 pattern</i> Colors <input type="checkbox"/> Berry <input type="checkbox"/> Brown <input type="checkbox"/> Gray	mediven 550 Design Elements <input type="checkbox"/> Animal <input type="checkbox"/> Crosses <input type="checkbox"/> Ornaments <input type="checkbox"/> Stars <input type="checkbox"/> Ribs <input type="checkbox"/> Pyramids				
mediven 550 Fashion Elements* <i>Must pick 1 color & 1 pattern</i> Colors <input type="checkbox"/> Berry <input type="checkbox"/> Brown <input type="checkbox"/> Gray	mediven 550 Design Elements <input type="checkbox"/> Animal <input type="checkbox"/> Crosses <input type="checkbox"/> Ornaments <input type="checkbox"/> Stars <input type="checkbox"/> Ribs <input type="checkbox"/> Pyramids							
<table border="1" style="width: 100%; border-collapse: collapse;"> <tr> <td style="width: 50%;"> Compression <input type="checkbox"/> CCL I (18-21 mmHg) <input type="checkbox"/> CCL II (23-32 mmHg) <input type="checkbox"/> CCL III (34-46 mmHg) (550 only) </td> <td style="width: 50%;"> <input type="checkbox"/> Hallux ease (550 only; seamless toe cap only) </td> </tr> </table>			Compression <input type="checkbox"/> CCL I (18-21 mmHg) <input type="checkbox"/> CCL II (23-32 mmHg) <input type="checkbox"/> CCL III (34-46 mmHg) (550 only)	<input type="checkbox"/> Hallux ease (550 only; seamless toe cap only)				
Compression <input type="checkbox"/> CCL I (18-21 mmHg) <input type="checkbox"/> CCL II (23-32 mmHg) <input type="checkbox"/> CCL III (34-46 mmHg) (550 only)	<input type="checkbox"/> Hallux ease (550 only; seamless toe cap only)							
<table border="1" style="width: 100%; border-collapse: collapse;"> <tr> <td style="width: 50%;"> Toe Options <input type="checkbox"/> Open Toes <input type="checkbox"/> Closed Toes <input type="checkbox"/> without Small Toe </td> <td style="width: 50%;"> Lymphpad (Optional) _____ Width (cm) _____ Length (cm) _____ Location </td> </tr> </table>			Toe Options <input type="checkbox"/> Open Toes <input type="checkbox"/> Closed Toes <input type="checkbox"/> without Small Toe	Lymphpad (Optional) _____ Width (cm) _____ Length (cm) _____ Location				
Toe Options <input type="checkbox"/> Open Toes <input type="checkbox"/> Closed Toes <input type="checkbox"/> without Small Toe	Lymphpad (Optional) _____ Width (cm) _____ Length (cm) _____ Location							

*Requires an extra 5 days for production.

medi USA

6481 Franz Warner Parkway
Whitsett, NC 27377

mediusa.com | phone 1.800.633.6334 | fax 1.888.570.4554