

JOBST® Custom Seamless Soft and Bellavar®



**Custom-made circular
knit for edema
management and
vascular conditions**

Lower extremity compression garments



Comfort, Health and Style!

JOBST® Custom Seamless Soft

JOBST® Custom Seamless Soft is the premium custom-made compression garment for mild to moderate vascular disorders. These garments are manufactured using precision technology and exact patient measurements. The use of attractive, lightweight fabric to create an anatomically fitting compression garment provides effective management of vascular indications and a stylish look.

Features and benefits:

Effective management for a wide range of indications

- Available in CCL 1 (18-21 mmHg*), CCL 2 (23-32 mmHg*) and CCL 3 (34-46 mmHg*)

JOBST Advanced Comfort®

- Added softness, moisture wicking to keep skin cool & dry

Multi-fiber double covered inlay yarns

- Increased softness & improved donning

Improved fabric distribution

- Enhanced fit and comfort with smooth appearance

Refined opacity

- Covers blemishes and scars

High stretch core yarns

- Increased recoil (stretch & recovery) to prevent sliding

Reinforced foot zones

- Increased strength & durability with added comfort

Optimized knee band

- Improved broader knitting construction and welt design; for added comfort and reduced risk of sliding / biting


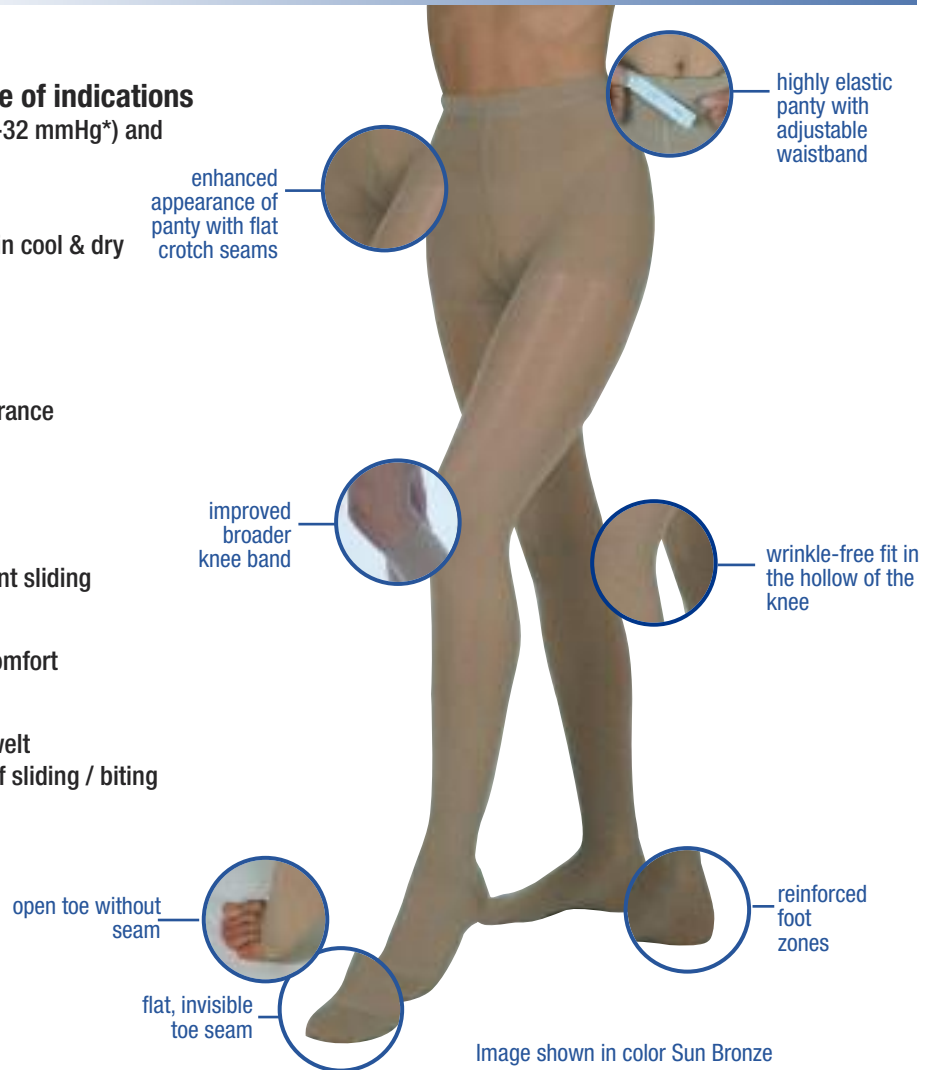


Image shown in color Sun Bronze

Improved Thigh Band Options



Soft knitted band (no silicone)

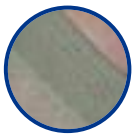

Silicone lace band



Silicone soft band (striped)



Silicone dotted band

Color Options




* Design Pressure
Federal (U.S.A.) Law restricts this device to sale by or on the order of a physician.

JOBST® Bellavar®

JOBST® Custom Bellavar® is the firm Custom garment for moderate to severe swelling of lower extremities. Bellavar is engineered with firm fabric to provide controlled stretch and resistance for therapeutic effect.

Features and benefits:

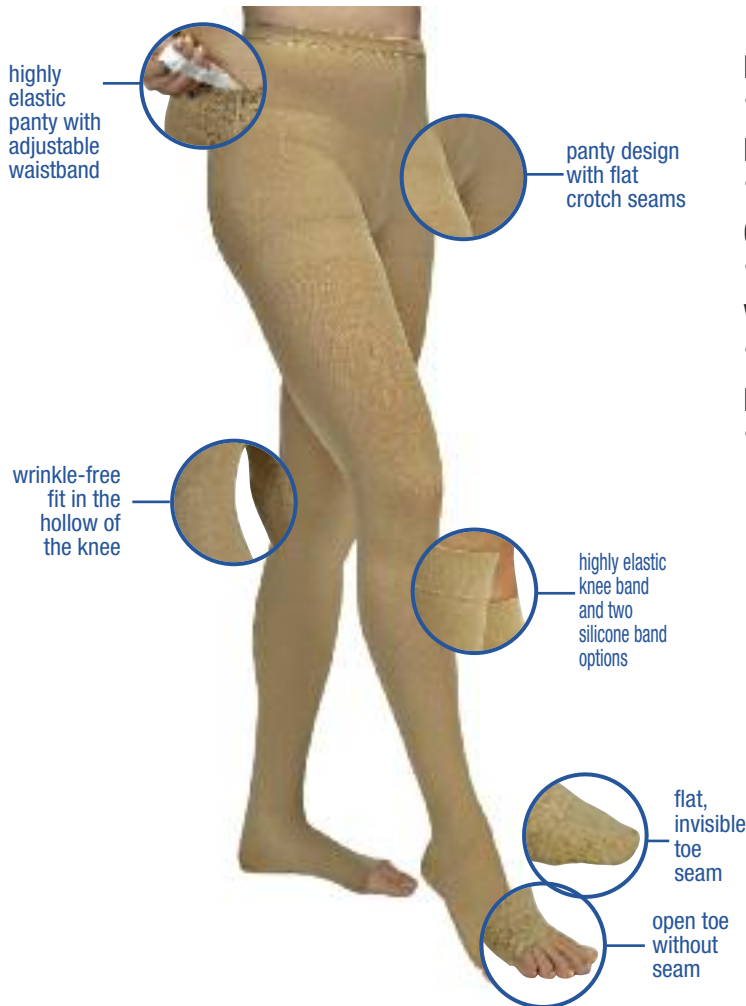

Firm fabric

- Provides controlled stretch and resistance for therapeutic effect

Effective management of a wide range of indications

- Available in CCL 2 (23-32 mmHg*) and CCL 3 (34-46 mmHg*)

Circular knit technology

- Creates a seamless garment

Variety of styles and options to satisfy patient needs

- 3 colors and 2 thigh band options

Easy measuring system

- Requires 30% fewer measurements than a conventional seamed garment

Basic Styles:

Knee Style



Thigh Style



Chaps Style



Thigh Band Options



plain thigh band with silicone



dotted silicone for plain thigh band



lace thigh band with silicone



silicone strips for lace thigh band

Color Options



Sand



Sun Bronze



Black

JOBST® Seamless Soft and Bellavar

Indications:

Compression Class I:

Mild swelling, tired aching legs, thrombosis prophylaxis, minor-varices, post-surgery

Compression Class II:

Primary varices, varices with edema, varices with pregnancy, post-sclerotherapy, post-venous ulceration

Compression Class III:

Secondary varices, severe varices with edema, post-phlebotic syndrome, chronic venous insufficiency, phlebitis; thrombophlebitis and edema

Contraindications:

- Dermatitis in the acute phase
- Arterial insufficiency
- Ischemia (e.g. advanced arterial disease)
- Wet dermatosis
- Uncontrolled congestive heart failure
- Untreated septic phlebitis
- Plegmasia cerulea dolens
- Cutaneous infections

Cautions / Recommendations:

- Incompatibility to fabric
- Impaired sensitivity of the limb
- Immobility (confinement to bed)
- Open wounds should be covered with an appropriate dressing before the compression garment is put on
- Ill-fitting compression garments may cause skin necrosis, pressure related peripheral nerve damage and deep vein thrombosis
- Conditions in which increased venous and lymphatic return is not desired

JOBST® Custom Seamless Soft and Bellavar®

Custom-made circular knit lower extremity compression garments

Ordering information:

Product No.		Description	Style
Seamless Soft	Bellavar®		
711020	711021	Knee High, Open Toe	
711040	711041	Knee High, Closed Toe	
713020	713021	Thigh High, Open Toe	
713040	713041	Thigh High, Closed Toe	
714020	714021	Chaps Style, One Leg, Open Toe	
714040	714041	Chaps Style, One Leg, Closed Toe	
715320		Chaps Style, Two Legs, Open Toe	
715340		Chaps Style, Two Legs, Closed Toe	
715520	715521	Waist High, Two Legs, Open Toe	
715540	715541	Waist High, Two Legs, Closed Toe	
716320	716322	Waist High, With Fly, Open Toe	
716340	716342	Waist High, With Fly, Closed Toe	
717520	717521	Maternity, Two Legs, Open Toe	
717540	717541	Maternity, Two Legs, Closed Toe	
718520	718523	Waist High, Full Compression, Open Toe	
718540	718543	Waist High, Full Compression, Closed Toe	
719520	719521	Waist High, One Leg Panty, Open Toe	
719540	719541	Waist High, One Leg Panty, Closed Toe	
	716321	Waist High, One Leg Panty, With Fly, Open Toe	
	716341	Waist High, One Leg Panty, With Fly, Closed Toe	
899115	899115	Color – Sand	
899111	899111	Color – Sun Bronze	
899116		Color – Amber†	
899113		Color – Navy†	
899114	899114	Color – Black	
799102	799102	5 cm plain/dotted silicone band (thigh / knee)	
799103	799103	6 cm lace / strips silicone band (thigh)	
799104	799104	6 soft / strips silicone band (thigh)	
799105	799105	2.5 cm plain / dotted silicone band (knee)	
799107	799107	Open pubis (waist high)	
799109		5 cm waistband (panty†)	
799115		2.5 cm waistband (panty†)	

† Not available in Bellavar

These products are measured using the landmark measuring process.

Please refer to the product label and package insert for full instructions on the safe use of these products.



BSN medical Inc. 5825 Carnegie Blvd. Charlotte, NC 28209-4633
 Tel. (+1) 704 554 9933 Fax (+1) 800 835 4325 www.jobst.com
 To order toll-free: JOBST (+1) 800 537 1063

60583 R1

© 2013 BSN medical Inc.

physician recommended



Comfort, Health and Style!

REV 02/13