

TributeNight[™] **Custom Lymphedema Garments** for Evening & Nighttime Therapy



www.Lohmann-Rauscher.us



In the early 2000's, lymphedema treatment was still being developed in the United States. The "gold standard" of treatment — multi-layer bandaging led to successful patient outcomes. However, patients often struggled to maintain therapy success on their own using bandaging. Lymphedema Therapists needed home therapy solutions for their patients managing this lifelong condition.

The Solaris Collection[™] answers therapists' challenges by providing effective solutions to help patients live with lymphedema and edema. It not only supports therapists and fitters, but helps improve patients' quality of life.

The Solaris Collection's range of edema and lymphedema treatment solutions will continue to expand, by listening to your clinical challenges and creative ideas, as well as your patients' needs.

- The L&R USA Family

Contents

About TributeNight [™]		
Powered by Coolcore®	6	
SURE-FIT [™] Guarantee		
Customer Support	10	
Build a TributeNight	13	
Style	14	
Channeling	15	
Profile	16	
Color	17	
Modifications	18	
Accessories	20	
How to Order		
Meet Our Models	24	
Indications		
Care Instructions		

About L&R USA INC.

L&R USA is an affiliate of our global parent company Lohmann & Rauscher, which has a 160-year history as a well-recognized and respected global manufacturer and supplier of compression therapy solutions and innovative wound care products. L&R USA provides you with effective compression solutions, including bandaging and the Solaris Collection[™] of garments and wraps, to help you treat your lymphedema and edema patients.

When you use L&R USA products, you know your patients are in good hands with trusted brands like TributeNight[™], ReadyWrap[™], Rosidal[®] K, and Debrisoft[®].

L&R USA — passionate about delivering innovative wound care and compression solutions.



About TributeNight[™] Your solution to maintain therapy success.

TributeNight[™] custom therapeutic garments provide lymphedema and edema therapy to patients during the evening and night. Designed by lymphedema therapists, TributeNight garments replace nightly multi-layer bandaging and offer an easy solution to maintain success.

TributeNight addresses your clinical challenges by offering a solution to meet your patients' needs to help improve their quality of life with lymphedema.

Therapy Challenges	TributeNight [™] Solutions
One size garment does not fit all •••••••>	Customized therapy for your patients' unique presentations
Reduced volume is not maintained> during the evening and night	Provides maintenance therapy while your patients are at rest and can help further reduce volume
Bandaging is difficult to apply correctly $\cdots $	Easy to don and doff appropriately
Lack of compliance with nightly bandaging \dots	Increases compliance by helping save time and reduce stress
Fibrotic tissue is a common issue that	Foam technology helps break up hardened, fibrotic tissue
Lack of compliance because	Powered by Coolcore [®] fabric to rapidly wick away sweat and moisture to regulate body temperature
Multi-layer bandaging does not •••••••••••••••••••••••••••••••••••	Allows freedom to move from evening throughout the night and during low activity
Family or caregivers are often>	Designed specifically for patients to apply easily on their own
Rolling bandages after washing	Easy garment care in the washing machine and dryer



Rest assured. TributeNight provides comfortable therapy through the night, while allowing your patients to rest.


TributeNight shown in Slate (Measurements from A–D)

You spoke. We listened.

Patients and therapists alike have asked for an evening and nighttime garment that offers a cool and comfortable experience - that is also clinically effective. That's why all TributeNight garments are now powered by Coolcore[®]!

Coolcore[®] is an innovative fabric with chemical-free Cooling Power which rapidly wicks away sweat and moisture to regulate body temperature. While other brands apply chemicals to the material for its cooling properties, Coolcore technology is built into the fabric to provide chemical-free cooling that makes it safer for your patients and the environment. This chemical-free technology maintains a lasting cooling benefit that withstands wash and wear.

Additionally, Coolcore fabric is soft, skinfriendly, and resists pills and snags, which offers your patients a comfortable garment experience for the one-year life of the garment.

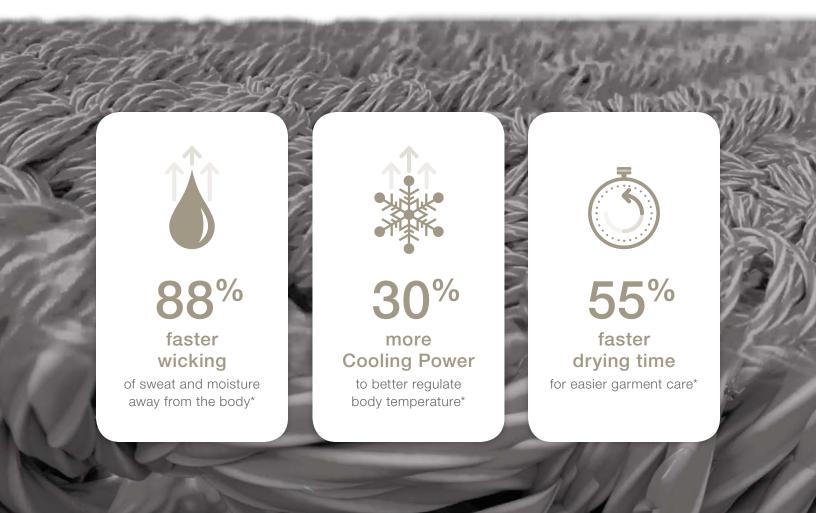





Moisture wicking

Sweat and moisture is rapidly wicked away from the skin, helping your patients stay dry.



Regulated evaporation

Patented fabric technology regulates the evaporation of moisture.



Moisture transportation

Moisture absorbed into the fabric is distributed across the surface of the fabric.



Chemical-free cooling

Cooling technology is built into the fibers of the fabric making it safer for your patients and the environment.



Resists pills and snags

Maintains garment integrity over its one-year life.



SURE-FIT Guarantee[™] Be confident your patient will have the right fit. Guaranteed.

We understand that measuring for custom garments is a combination of art and science. As part of our promise of service excellence to you, we want to ensure you are confident in the design and customized fit of each TributeNight.

To support you and your patients, all TributeNight garments are covered, at no extra charge, by our SURE-FIT Guarantee[™]. This complimentary guarantee is comprised of three industryleading services for the full life of each TributeNight garment: First Fit, Reduction Alteration, and Materials & Manufacturing.



1

First Fit

1-45 DAYS

To ensure the proper fit, we will alter or re-make your patient's garment as many times as needed within the first 45 days from shipment.

2

Reduction Alteration

45 DAYS-1 YEAR

To support your patient's successful compliance, we offer one free alteration during the one-year life of the garment to reduce the size of the garment to match their improved measurements.

3

Materials & Manufacturing

1 YEAR

We guarantee your high-quality garment to be free of material and manufacturing defects during the recommended one-year life of the garment.

Customer Support Supporting you and your patients, every step of the way.

Custom Design Center

We offer therapists and fitters direct access to our team of designers, to lend their expertise in helping design custom garments for patients from coast to coast by:

- answering measuring and ordering questions
- collaborating with you on garment design for your patient's unique presentation and abilities
- reviewing patient order history and current measurements for accurate reordering
- redeeming SURE-FIT Guarantee
- making recommendations for fit, style, and modifications

PASSION IN THE MAKING

43

Combined years of TributeNight custom design and manufacturing support for therapists and fitters

Years designing and manufacturing TributeNight garments in Milwaukee, WI USA

Guaranteed business days from receipt of order to arrival at requested delivery address*

4

The Custom Design Center is available Monday through Friday, 7:00AM-7:00PM Central Time

- **** 414-892-5158
- CustomDesignCenter@us.LRmed.com
- Video conference
- Travel







Patient Outreach Service

L&R USA INC. extends service support directly to your patients through our Patient Outreach Service. We are available to address all of your patients' garment needs, including:

- how to register their garment to activate the SURE-FIT Guarantee
- how to properly don and doff
- garment care

Available Monday through Friday, 7:00AM–7:00PM Central Time

& 414-892-4140

PatientOutreach@us.LRmed.com

TributeNight shown in Purple (Measurements from A–G) Alti i i

Build a TributeNight. Completely customizable.

Every lymphedema case expresses uniquely, from the affected area to the various presentations in shape and size.

Because of this, the best maintenance garments are customized to your patients' exact measurements to provide the right support for your patients' therapy needs.

TributeNight garments are custom designed and individually made for your patients, incorporating their style, fit, functionality, and donning preferences.

The following pages highlight the range of custom design options, but does not represent every way TributeNight can be customized.

Both you and your patient deserve a garment to support your success.

Building a TributeNight is easy.

Follow these simple steps.



Explore the customization possibilities with TributeNight >



Determine garment style. Coverage from head to toe and in-between.

TributeNight is customizable to address your patients' lymphedema wherever it presents itself. To begin, identify where you would like the garment to cover, then use our measuring form to indicate the length and circumference measurements where coverage is needed.

Below are some ideas to get you started:






Select the channeling type. Designed based on MLD concepts.


TributeNight garments have deep channels sewn into the garment, which can be ordered in a diagonal chevron pattern or straight vertical pattern.

These channel designs are based on the Manual Lymphatic Drainage (MLD) therapy concepts of helping direct the flow of built-up lymphatic fluid.

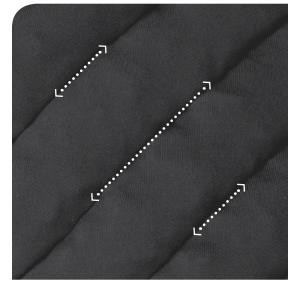

Chevron channeling

- Designed to direct lymphatic fluid to collateral pathways.
- For patients with secondary lymphedema.

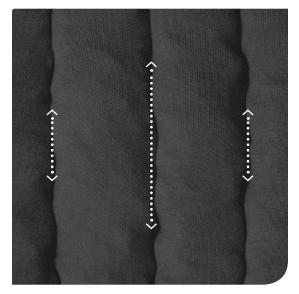


Vertical channeling

- Designed to direct fluid to intact lymph nodes and the primary lymphatic system.
- For patients with primary lymphedema.




Choose profile height. Offering foam technology variety benefiting patient comfort.

TributeNight's foam technology helps break up fibrotic tissue to assist in the softening of skin to help improve lymphatic flow. The TributeNight has a blend of high and low density chopped foam that creates multi-point pressures on the affected area covered by the garment. The foam technology is offered in two different profile heights based on you and your patient's preference – Original and Low profile.

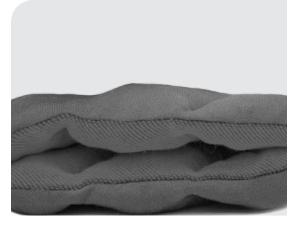

Original profile

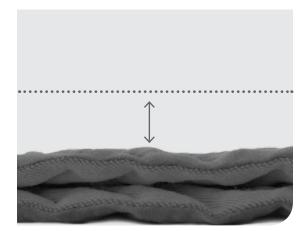

Low profile





Pick fabric color. Patients can express themselves.

Our range of fashionable garment fabric colors offers your patients the opportunity to participate in designing their garment.





Add optional modifications. Optional features to enhance patient compliance.

Patients with limited mobility and dexterity benefit from modifications which tailor the garment to their abilities. Some modifications assist with donning and doffing such as zippers and pull-up loops, while others promote safety, like our non-skid pads.

These modifications can be added to any TributeNight design.

- Zippers
- Closure, 10cm maximum (VELCRO® brand)
- Adjustable panels (VELCRO® brand)
- Pull-up loops
- Snap tape
- Digit spacers
- Non-skid pads
- Lip bridge
- Tracheotomy accommodations







Digit Spacers









Add accessories. Enhance the garment and their experience.

TributeNight garments can be paired with optional accessories to enhance your patients' garment experience. These accessories are designed to provide better containment, protect garment from external elements or simply support easy donning.

Variable Compression Jacket (VCJ)

Helps provide extra support and containment for patients with hard-to-manage and highlyfluctuating lymphedema and edema.

Adjustable VELCRO® Brand fasteners can be preset, allowing patients to don over their TributeNight and secure with an easy-to-pull zipper.



Outer Jacket (OJ)

Helps provide support and containment, while protecting the garment from dirt, pet hair, and stains.

Choose VELCRO® Brand or snap fastener. (See pages 18 and 19)




Donning Aid

Easy-Slide

Easy-Slide donning aid is available to support selfcare by providing an easy don & doff experience.



How to order TributeNight. Easy steps for custom therapy.

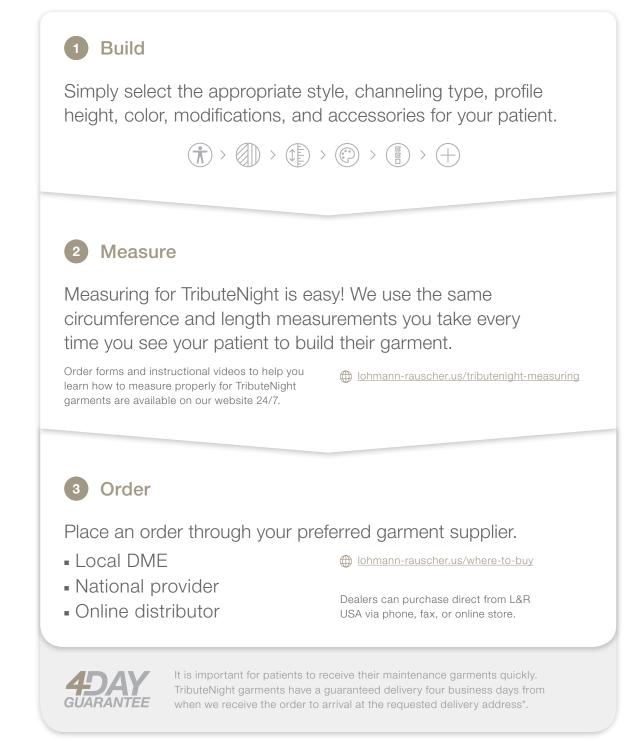

*Orders must be received before 2:00 PM Central, Business days only, Holidays excluded, US customers only, On hold garments receive guarantee starting at release of the on hold garment, Offer available for new orders only (warranty garments do not apply), Offer available for all TributeNight custom, ExoCustom, and ReadyWrap custom orders.



Have questions? Our customer support team is here to help.

Custom Design Center

The Custom Design Center is available Monday through Friday, 7:00AM–7:00PM Central Time

& 414-892-5158

⊠ CustomDesignCenter@us.LRmed.com

□ Video conference

⑦ Travel

Patient Outreach Service

Our Patient Outreach team is available Monday through Friday, 7:00AM–7:00PM Central Time

& 414-892-4140

PatientOutreach@us.LRmed.com



Meet our models: patients living with lymphedema.

Becca, 36

"I love TributeNight. It is custom to you, because you are different than everyone else. They provide me comfortable therapy at night, are easy to don on and off, durable, and allow great movement. My nightly care for my Lymphedema is less time consuming. I would never go a night without my TributeNight garments, they are a part of my everyday life."

ABOUT

Becca is a wife, and mother of five children

HOBBIES

Co-Chair of the Wisconsin Chapter of the Lymphatic Education & Research Network (LE&RN), active in the Lymphedema Treatment Act, oil painting, crocheting, anything crafty or artistic

LYMPHEDEMA

Primary in both legs, both arms, back, and trunk


DIAGNOSED IN 1999, after nine years undiagnosed

DAYTIME GARMENTS

ExoCustom[™] thigh highs and knee highs, ExoCustom arm sleeves, and ExoStrong[™] gauntlets



TributeNight shown in Purple (Measurements from A–G)

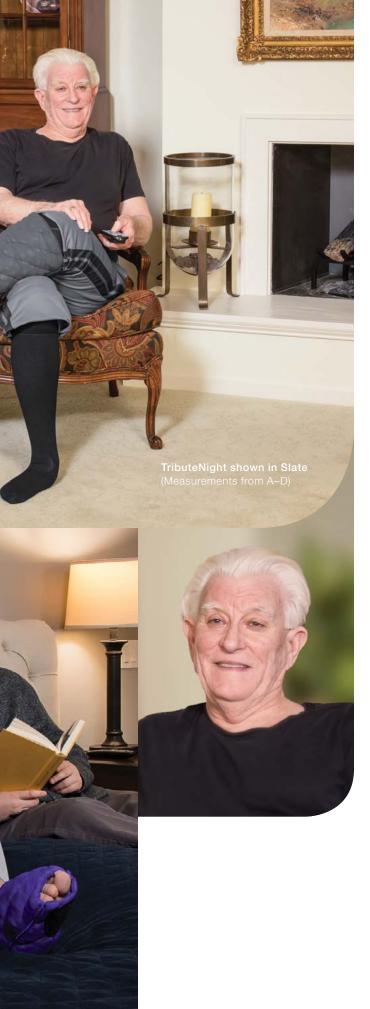

Peter, 74

"TributeNight is the best night garment I've ever worn, because it is comfortable and breathable. It also goes on so easy. You wake up in the morning, and when you get out of bed your swelling is maintained."

ABOUT PETER

Father of four, grandfather of five, Vietnam War Veteran

HOBBIES Fishing, working out, reading

LYMPHEDEMA

Secondary in lower extremity after a surgery

DIAGNOSED IN 2015, previously misdiagnosed with edema

DAYTIME GARMENT

ReadyWrap™

"The ReadyWrap helped me turn a corner with my lymphedema. I wear ReadyWrap when I get up in the morning until I get ready for bed with my TributeNight."

Indications

INDICATIONS

- Various stages of lymphedema
- Edema (non-blood related)

INTENDED USES

- Helps maintain reduced volume and shape by supporting extremities and torso
- Management of various stages of lymphedema
- Management of edema (non-blood related)
- Assists in softening of indurated/fibrotic tissue

CONTRAINDICATIONS

- Arterial insufficiency or ABPI < 0.8
- Moderate to severe peripheral arterial disease
- Untreated congestive heart failure
- Untreated infection (e.g., cellulitis)
- Allergies or sensitivities to materials

PRECAUTIONS

- Impaired sensation (peripheral neuropathy)
- Cognitive impairment

Care Instructions

WASHING

- Machine wash on permanent press in warm water with other laundry items
- Do not use bleach or fabric softener
- Take extra care if the garment has VELCRO® Brand fasteners and ensure it is fastened securely prior to washing
- Do not dry clean

DRYING

- Tumble dry on permanent press with other laundry items or clean tennis balls
- Drying time may be shortened by pressing out excess water with a towel prior to drying
- Do not iron



MATERIAL CONTENT

Polyurethane, polyester, spandex

Accessory information available upon request.

TributeNight shown in Purple (Measurements from A–G)





L&R USA INC. 3880 W Wheelhouse Road Milwaukee, WI 53208 USA
 Phone:
 855-892-4140, 414-892-4140

 Fax:
 414-892-4150

 Email:
 inquiries@us.LRmed.com

Debrisoft[®] and Rosidal[®] are registered trademarks of Lohmann & Rauscher. Coolcore[®] is a registered trademark of Coolcore, LLC. VELCRO® is a registered trademark of Velcro BVBA.